

EVIS EXERA III

EVIS EXERA III Bronchoscope

$BF-1TH190 \\ \label{eq:Fully rotatable}$ Fully rotatable, therapeutic bronchoscope with superb HDTV image quality.



EVIS EXERA III Bronchoscope

BF-1TH190



The revolutionary design of the BF-1TH190 bronchoscope, powered by EVIS EXERA III, advances therapeutic bronchoscopy with outstanding handling, improved functionality, and unmatched image quality.

Advancing Maneuverability

With outstanding handling and unmatched tracheobronchial access

- Rotary Function: Insertion tube rotation of up to 120° left and right may effectively reduce hand rotation* and supports easier insertion and positioning of endoscopic devices.
- Waterproof One-touch Connector: A simple one-touch connector minimizes the effort required for set up, and its fully submersible design eliminates the need for a water-resistant cap.




*When compared to the BF-180 series.

Advancing Visualization

With HDTV for clearer, brighter imaging

- HDTV: The brighter, crystal clear images delivered by true HDTV enable more detailed observation of bronchial surfaces during bronchoscopy.
- Electronic Magnification: Close-up observation of bronchial surfaces is readily available with 1.2 x and 1.5 x electronic magnification.
- Narrow Band Imaging (NBI): Now significantly brighter compared to previous generation scopes, NBI provides twice the viewable distance for enhanced visualization of vascular structures and mucosal morphology.

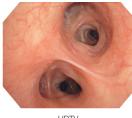

HDTV



NBI

Advancing Versatility

With compatibility across multiple specialties

- EVIS EXERA III: The power of EVIS EXERA III expands across the hospital to enable equipment standardization and simplified support and service, with the ultimate goal of lowering the total cost of equipment ownership.
- Endoscopic Devices: A full array of Olympus bronchoscopy devices are available and designed to deliver superior performance while minimizing channel damage, thereby potentially reducing repair costs and extending endoscope life.

Specifications		
Optical System	Field of View	120°
	Direction of View	Forward Viewing
	Depth of Field	3 to 100 mm
	Image Quality	HDTV
Insertion Tube	Distal End Outer Diameter	6.2 mm
	Distal End Enlarged	
	Instrument Channel Outlet	Light-guide Lens
	Right Left Objective Lens	
	Insertion Tube Outer Diameter	6.0 mm
	Working Length	600 mm
	Rotary Function	Yes
Instrument Channel	Channel Inner Diameter	2.8 mm
	Minimum Visible Length	3.0 mm from the distal end
	Direction from which endoscopic devices enter and exit the endoscopic image	
Bending Section	Angulation	Up 180°
Total Length	880 mm	
Electrocautery Instrument Compatibility	Yes	
Laser Compatibility	Nd: YAG, 810 NM diode	
Compatible EVIS EXERA System	Video System Center Xenon Light Source	Olympus CV-190 Olympus CLV-190

Olympus is a registered trademark of Olympus Corporation, Olympus America Inc., and/or their affiliates. | Medical devices listed may not be available for sale in all countries.



OLYMPUS AMERICA INC. 3500 Corporate Parkway, PO Box 610, Center Valley, PA 18034 For more information, contact your Olympus sales representative, or call 800-848-9024. www.medical.olympusamerica.com

©2018 Olympus America Inc. All rights reserved. Printed in the USA OAIRES1018BRO28704