

OLYMPUS[®]

Your Vision, Our Future



EVIS EXERA GASTROINTESTINAL VIDEOSCOPE
OLYMPUS GIF TYPE 1TQ160

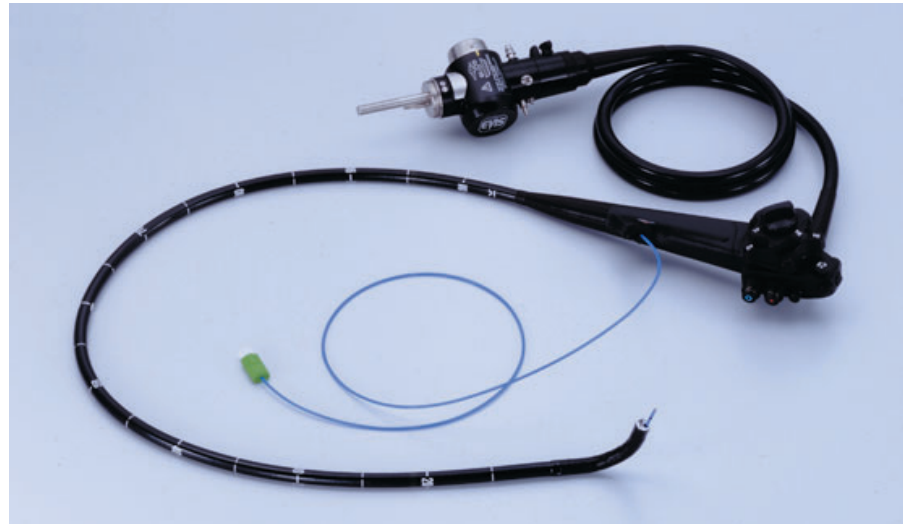
LARGE CHANNEL AND EXCELLENT IMAGE QUALITY

Take advantage of large-screen,
high quality image in a versatile videoscope with
a large 3.7 mm channel design
for extended instrumentation capability
even while suctioning




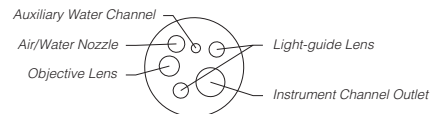
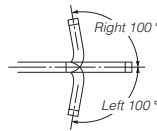
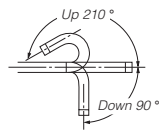
Main Features for Versatile Performance

- The large 3.7 mm diameter channel incorporated in the slim 11.3 mm insertion tube not only enables use of a wide range of instruments, it even allows suctioning with an instrument inserted.
- Large display size delivers outstanding quality with sharper, crisper images and life-like color from edge to edge.
- Extra-wide 140° field of view enables accurate observation of a wider area.
- When used in combination with our dedicated flushing pump and CV-160 video system center, a continuously clear view is maintained by a forward water jet that removes mucus and clots inside the upper digestive tract at the touch of a scope switch. The connection for this function has been moved to the light-guide connector section to keep it out of the way of the operator's hand.
- Four-way angulation (210° up, 90° down and 100° right/left) enables comprehensive examination of the upper digestive tract.
- Ergonomically designed grip enhances scope maneuverability while redesigned, easy-to-access knobs and user-selectable switches improve operability.
- Fully compatible with the CV-160/140.
- Scope ID function stores individual scope information and displays it on the monitor to facilitate endoscopy suite management.



Specifications

Optical System	Field of view	140°
	Direction of view	Forward viewing
	Depth of field	3 to 100 mm
Distal End	Outer diameter	10.9 mm
Insertion Tube	Outer diameter	11.3 mm
Bending Section	Angulation range	Up 210°, Down 90°, Right 100°, Left 100°
Working Length		1030 mm
Total Length		1340 mm
Instrument Channel	Inner diameter	3.7 mm
	Minimum visible distance	4 mm from the distal end
	Endotherapy accessory entrance/exit position in field of view 	



Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.

OLYMPUS[®]

OLYMPUS MEDICAL SYSTEMS CORP.
Shinjuku Mondlioth, 3-1 Nishi-Shinjuku 2-chome, Shinjuku-ku, Tokyo 163-0914, Japan
OLYMPUS MEDICAL SYSTEMS EUROPA GMBH
Postfach 10 49 08, 20034 Hamburg / Wendenstrasse 14-18, 20097 Hamburg, Germany
OLYMPUS AMERICA INC.
3500 Corporate Parkway, P.O. Box 610, Center Valley, PA 18034-0610, U.S.A.
OLYMPUS LATIN AMERICA, INC.
5301 Blue Lagoon Drive, Suite 290 Miami, Florida 33126-2097, U.S.A.
KEYMED LTD.
Keymed House, Stock Road, Southend-on-Sea, Essex SS2 5QH, England

OLYMPUS SINGAPORE PTE LTD.
491B River Valley Road #12-01/04, Valley Point Office Tower, Singapore 248373
OLYMPUS HONG KONG AND CHINA LIMITED
Room 1520-1527, Ocean Centre, 5 Canton Road, Tsimshatsui, Kowloon, Hong Kong
OLYMPUS (BEIJING) SALES & SERVICE CO., LTD.
Room.1202, NCI Tower, A12 Jianguomenwai Avenue Chaoyang District Beijing, 100022, China
OLYMPUS MOSCOW LIMITED LIABILITY COMPANY
117071, Moscow, Malaya Kaluzhskaya 19. bld. 1, fl.2, Russia
OLYMPUS AUSTRALIA PTY LTD
31 Gilby Road, Mount Waverley, VIC., 3149, Australia